





Chair Topper® eBrochure

Developed to assist wheelchair users who are able to transfer, the **Chair Topper** automatically folds and stores a conventional folding wheelchair inside a stylish, weathertight car-top compartment. The need to store your chair in the trunk or back seat is eliminated, allowing family and friends to enjoy the comfort of your vehicle. The universal mounting system ensures that it will fit most automobiles and light trucks. It adds only 55 cm in height, and the sleek plastic shell gives a finished, automotive look. A convenient hand-held control operates the folding, lifting and storing of your wheelchair. It takes only 30 seconds to complete the cycle, and the smooth operation ensures your chair is carefully handled and transported without undue wear and tear.

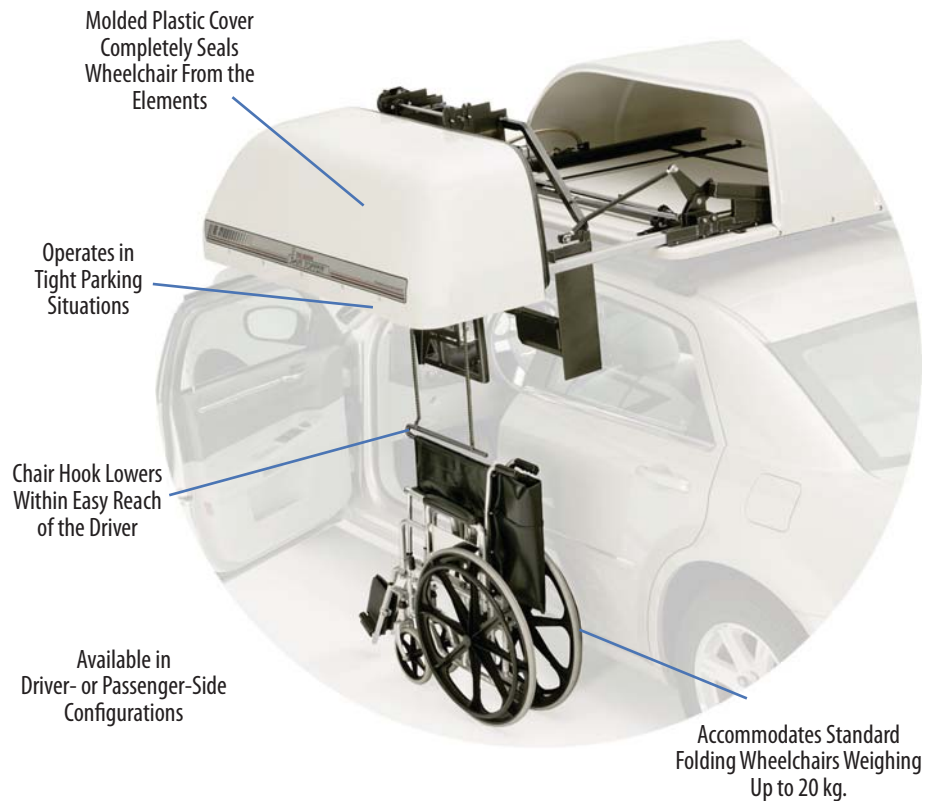

Car-Top Carrier for Folding Wheelchairs

Please note: this installation is shown as an example only. Mounting system shown is no longer available. Before purchasing a Chair Topper, please consult with your local BraunAbility Dealer to ensure it will fit on your vehicle.

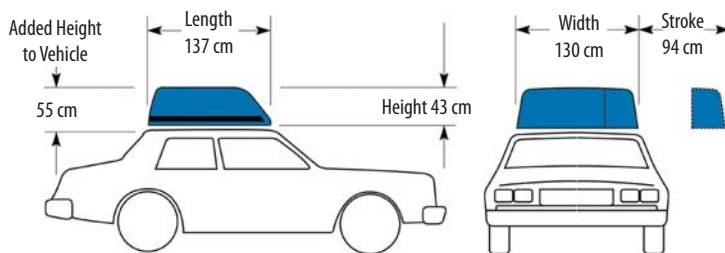

Please consult with your local BraunAbility Dealer for a list of vehicles that can accept the Chair Topper.



Electromechanical, plastic-shell car-top wheelchair carrier designed to fit a wide variety of domestic and foreign vehicles. Lifting Capacity is 20 kg. The Chair Topper utilizes a 12 volt DC reversible automotive-type motor, powered by the vehicle battery. The aluminum frame is encased in a molded plastic shell, which completely seals the wheelchair from the elements. The unit weighs 56 kg. The Chair Topper is shipped with a permanent mounting system. Contact your area sales manager to ensure it will fit on your vehicle. Also, please refer to the specification information in this guide before purchasing a Chair Topper. Special considerations may be required for vehicles with luggage racks, sun roofs, hatchbacks, etc.



+1 574 946 6153

www.BraunAbility.com/international

All illustrations, descriptions and specifications in this brochure are based on the latest product information at the time of publication. The Braun Corporation reserves the right to make changes at any time without notice. © 2010 The Braun Corporation